

## Optional accessories for individual requirements



The manually adjustable footrest is designed to support the legs in sitting position.



With additional handles the chair can be controlled from the foot end.




The dinner tray is clamped simply and securely onto the side rails.




One or two infusion poles can be mounted easily on receptacle systems.



The receptacle for the remote control allows clear and direct operation.

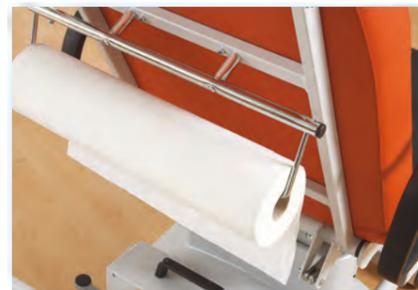

The paper roll holder allows for increased infection control and hygiene.

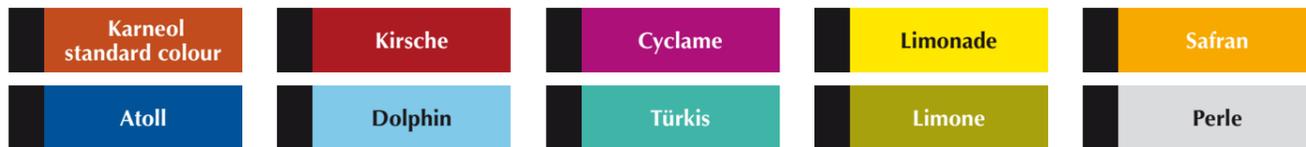

**Colouring:** The basic colour of the mattress is always black combined with a second colour. Standard version is black-kaarneol as seen on the pictures. There are many more colours available. Please ask for the reliable colour chart.

### Technical Data

#### Dimensions

Seat segment width	60 cm
Chair width incl. guard rails	78 cm
Total length	192 – 210 cm
Total weight	106 kg

#### Adjustment features

Back segment angle	80 – 0°
Seat segment angle	30 – 5°
Leg segment angle	50 – 0°
Height adjustment	63 – 88 cm

#### Electrical specifications

Protection category	IPX 4
In compliance with	EN 60601-1
Motors	4

#### Safety

Safe maximum load	180 kg
-------------------	--------



Bionic Medizintechnik GmbH ... a JMS-corporation · Max-Planck-Straße 21 · D-61381 Friedrichsdorf/Germany  
Telefon: +49(0)6172-75760 · Fax: +49(0)6172-757610 · info@bionic-jms.com · www.bionic-jms.com



**TekMed™**

1300 720 727

12 Bellevue Crescent Preston Victoria 3072 Australia  
T (03) 9487 9999 F (03) 9484 8792  
info@tekmed.com.au www.tekmed.com.au



**TekMed NZ™**

Ph: 0800 456 217

info@tekmed.co.nz

www.tekmed.co.nz

Right for technical modifications reserved B-MOBRE/03 2013-08

## Flexibility and comfort in your recovery room



 **RecoLine<sup>+</sup>**

  
www.therapychair.com

# RecoLine+ – flexibility and comfort for your patients and nursing staff in your recovery room



The electrical height adjustment allows smooth and convenient transfer between the operation table and the RecoLine+



Obtain perfect height and positioning for treatment with just one button



Individual adjustment of all segments of the chair allows optimal seating positions and faster recovery for the patient



The upright seating position provides secure and convenient entry and exit before and after the ambulant surgery

**The RecoLine+ therapy chair has been developed to improve patient care after minor operations taking into consideration the needs of both patients and nursing staff, giving more flexibility and comfort in the recovery room.**

Four motors offer a continuous adjustment of both lying and sitting positions as well as height adjustment by an easy to operate remote control.

Large wheels with a diameter of 12,5 cm and foldable handles on the back segment allow safe and fast transportation of your patients.

A quickly rechargeable, strong battery with 24 V and 7,2 Ah capacity provides off-line operation for several days.

In case of an emergency situation the shock position can be activated by a red foot pedal. This enables the medical team to start emergency measures immediately without pressing the emergency button by hand first.

The guard rails can be folded partly or totally to facilitate the entry into the chair for the patient.

Numerous optional accessories are available to fit the RecoLine+ to your facility's individual requirements. Available are a dinner tray, additional handles on the leg segment, receptacle for infusion pole, a paper roll holder and a receptacle for the remote control.



The comfortable pillow provides a pleasant reclining position. The pillow can be positioned depending on the patient's size.



The guard rails can be folded partially or completely to facilitate the entry and the transfer of the patient.



The 12.5 cm wheels are equipped with a steering lock for easy driving and central brake for secure positioning.



The battery is easily rechargeable via a charging cable and guarantees an off-line operation for several days.



The shock position can be easily activated by a red foot pedal. Thus, the nursing staff has the hands free for further measures.



The three most important positions are pre-programmed on the remote control. For a fine-tuning of positions there are additional keys.